

# Keles Face Mask

## Patient Instructions

The intraoral appliance is activated  $\frac{1}{4}$  turn (one activation) twice daily for the first 5 days. The Keles Face Mask should then be worn for at least 14 hours daily, including sleep, for the next six months.

### To remove the mask:


- 1) Disengage one side of the elastic band on the forehead rest from the vertical wire bar.
- 2) Hold the horizontal bar and push it back into the mouth until the bow is disengaged from the tube on the molar band, approx.  $\frac{1}{4}$  inch.
- 3) Lower the facebow and remove from mouth.

### To replace the mask:


- 1) Secure the Keles Face Mask into the tubes of the intraoral appliance
- 2) Connect the forehead unit to the hook on the vertical bars with the elastics.

### Silopad Maintenance:

Use mild face soap & luke warm water. Do not remove pads to wash or use any alcohol-based soap. Blot dry with a clean cloth.

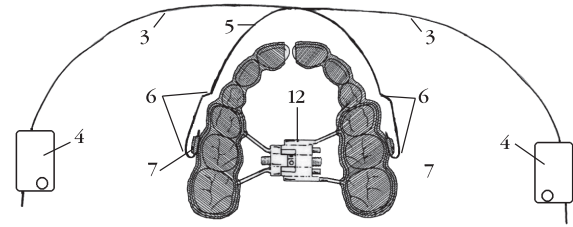

**Figure 1** Intraoral RPE appliance with Keles Face Mask. (3) Horizontal bar. (4) Adjustment blocks. (5) Intraoral bow. (6) Chairside bends to customize intraoral bow to tubes of RPE appliance. (7) Tubes to connect RPE appliance to Keles Face Mask. (12) RPE appliance.

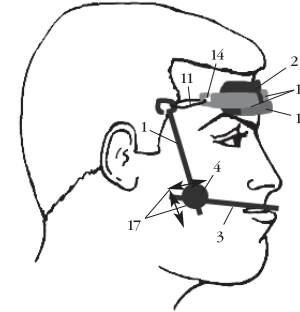

**Figure 2** Extraoral components of Keles Face Mask on patient. (1) Vertical bar with elastic hook. (2, 10, 13, 14) Adjustable forehead components. (3) Horizontal bar. (4, 17) Adjustment blocks with Allen screw locks and activation movement. (11) Elastic elements.