


The Original Adaptable Class III Mask

Standard and Pedo Size Available

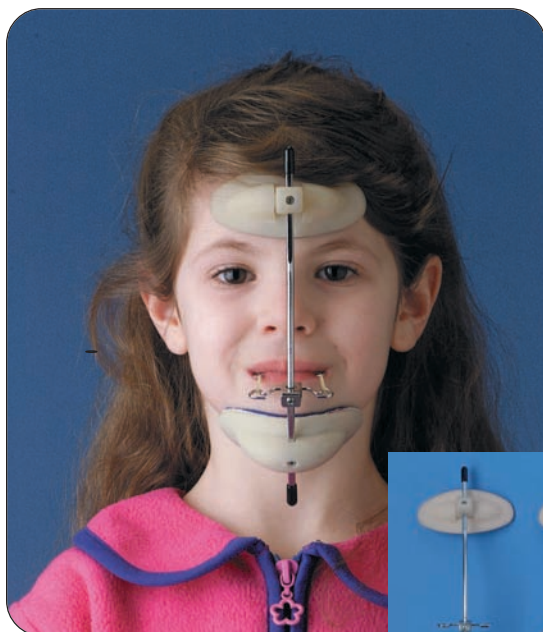

Order No.
HST - Standard Mask
HSM - Pedo Mask
Single use only

Adaptable Class III Mask

Dr. Henri Petit

- Made in U.S.A.
- Best Quality Materials and Workmanship
- Designed and Approved by Dr. Henri Petit

Gel liner for forehead and chin cup maximizes comfort (LATEX FREE)

Washable, reusable Silopad™ liner effectively reduces:

- Friction
- Abrasion
- Shear Force
- Pressure

Forehead Rest

- Patient comfort is maximized with ample coverage for stability and heavy force diffusion.
- Shatter-resistant high-impact ABS plastic.

Mouthbow

- Looped ends for patient safety.
- Radial locking block is repeatedly adjustable from low to high pull.
- Welded crossbar.

Chin Cup

- Peripheral design offers laterally extended coverage to diffuse force over the mandible while avoiding the soft area near the lip and under the chin.

Main Frame

- Arc of opening bend in lower portion simplifies adjustment.
- Length adjustable at both ends, utilizing the forehead rest Slide/Lok and chin cup guide position.
- Adjustable Pedo Size 7-1/2" in length.
- Adjustable Standard size 9-1/2" in length.

IMPORTANT REMINDER FOR DOCTORS!

To avoid risk of eye injury, be sure to inform patients and parents about proper safety precautions when placing and removing headgear. Also caution patients against rough or risky play while headgear is in place.

