

Great Lakes Operator Pliers



Distal End Cutter with Safety Hold- Tool Steel Insert

221-001

Cleanly cuts and securely holds wire. Angled tip allows arch wires to be cut close to the distal end of buccal tubes. Maximum cutting capacity: .0215" x .028" stainless steel. Smallest suggested wire size: .014" nickel titanium



Omega Tweed Loop Plier

221-008

Form uniform Omega loops in square and rectangular wire up to .0215" x .028." Each

section of the round beak has a length of .090" and diameter of .045", .060", and .075".

Replacement Omega Tips 221-009
(2 tips and wrench)

- Manufactured with high quality German stainless steel forgings to ensure consistent instrumentation
- Adjustable tension for greater handling control
- Superior life expectancy



Distal End Cutter

221-003

Cuts arch wire nearly flush with the distal end of the buccal tube,

helping to eliminate potential patient irritation. Excellent for wires smaller than .014" nickel titanium. Maximum cutting capacity: .0215" x .028" stainless steel



Tweed Rectangular Arch Plier

221-010

Smooth working edges will not nick or mar

archwire. Short jaws provide good control, and are parallel when opened to .020", ensuring accurate 90° bends.



Weingart Plier

221-012

Easy access to lingual and posterior areas is provided by tapered, slightly rounded beaks that fit easily between

brackets. Pointed, serrated pads are precision matched assuring better balance and firm grip.



Mini Pin and Ligature Cutter-Tool Steel Insert

221-004

Fine, durable tip for easy access to tight areas. Cuts the smallest soft wire, pins, ligatures, and elastics. Maximum cutting capacity: .015" soft wire



Sliding Hook Crimping Plier for Arch Hooks

221-011

Thinner beaks accommodate

curvature of arch hook. Firm grip on tube is maintained by two grooved inserts. The tube is securely crimped at three points to prevent the ends from flaring.



Posterior Band Remover-Long

221-013

Designed for easy access to distal areas. Perfect for posterior band removal when

padded tip is on occlusal of the molar cusp.

Plastic replacement pad,
3/16" diameter (1 pad) 221-014

Great Lakes Operatory Pliers Provide Optimal Performance!



Bird Beak Plier

221-015

For round wire up to .030", working edges are beveled to prevent nicking or marring of wire. Cone shaped beak opposed by "pyramid" shaped beak.



How Style Straight Utility Plier

221-020

Slimmer beaks are designed to grasp all types of auxiliary attachments. Precisely matched, serrated tips hold wire at any convenient working angle. Rounded beaks provide tissue safety and patient comfort.



Mathieu (Elastics)

221-022



Positive locking plier with extra-fine, serrated beaks for holding and placing ligatures and elastics without tearing.



Jarabak Plier

221-016

Precision grooves and serrated tips ensure uniform, accurate bending and closing of loops in lighter wire.



Three Prong Plier


221-021

Precision aligned tips bend consistently, and gradual tapering allows delicate bending. Single piece construction of the double tip eliminates flaring and ensures strength.



Mosquito Elastic Placing Forcep

221-024

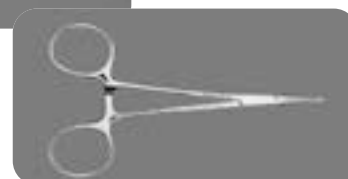

Very delicate, cross-serrated tips guaranteed not to tear elastics or allow them to slip.



Padded Direct Bond Bracket Remover-Wide

221-017

Pop-off brackets comfortably and easily.

Pad rests on incisal surface and chisel tip

wedges between bracket and base of tooth. Can also be used to remove remaining adhesive.

Replacement tool steel tips (2) 221-018
and wrench (1)

Plastic replacement pad, 221-014
3/16" diameter (1 pad)

To order call toll-free:

800-828-7626

Local: 716-871-1161

Fax: 716-871-0550

E-mail: info@greatlakesortho.com



Great Lakes

ORTHODONTICS, LTD.

An Employee Owned Company